


Pluggo RH Architect 5 Position Racks Automated Decapper - Technical Description

Pluggo is designed for the safe removal of original caps from most standard vacuum collection tubes. Pluggo can handle 5ml, 7ml and 10ml vacuum tubes with 75 x 13 mm, 100 x 13 mm, 100 x 16 mm tubes carried in the Architect position racks.

PLUGGO

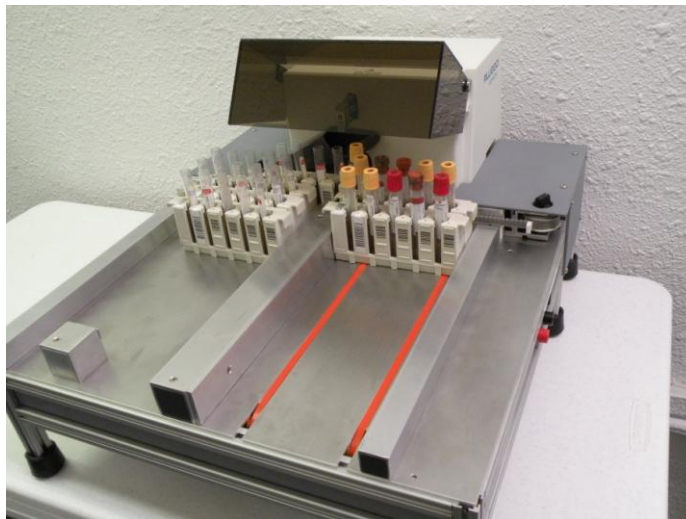

A detachable chute guides the caps into a disposal area beneath the instrument after decapping.

Technical data

Capacity:	Over 2,000 tubes per hour
Sample tubes:	Standard vacuum blood tubes
Size:	56 x 61 x 36cm (Depth x Width x Height)
Weight:	23 kg
Voltage:	230VAC 50Hz / 110VAC 60Hz
Power:	60VA