


## Pluggo Standard Automated Decapper - Technical Description

Pluggo is designed for the safe removal of original caps from most standard vacuum collection tubes. Pluggo is capable of handling 5ml, 7ml, 7ml long and 10ml vacuum tubes in the same machine, using different carousels.



# PLUGGO

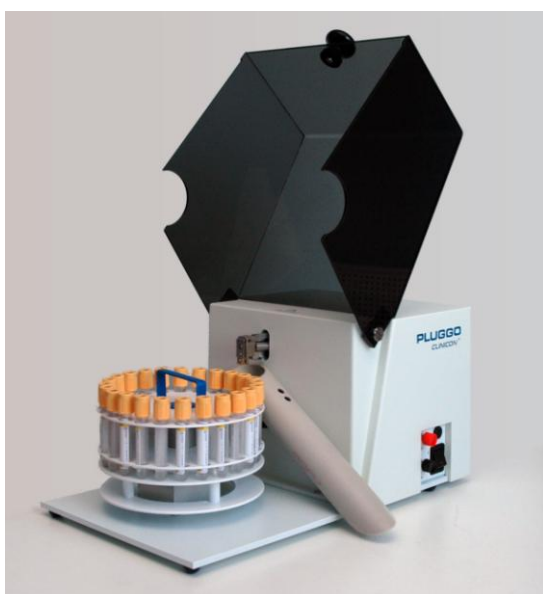

Carousel colour	Tube size	Volume
Blue	13 x 75 mm	5 ml
Red	13 x 100 mm	7 ml long
Grey	16 x 75 mm	7 ml short
Yellow	16 x 100 mm	10 ml

A detachable chute guides the caps into a disposal area after decapping.

### Technical data

Capacity:	Up to 24 samples in each carousel
Sample tubes:	Standard vacuum blood tubes
Speed:	30 seconds per carousel
Size:	43 x 27 x 21cm (Depth x Width x Height)
Weight:	9,0 kg
Voltage:	230VAC 50Hz / 110VAC 60Hz
Power:	60VA